

oriori®

Smart Grip Ball

motion series

Manual

POWERED BY
LIGHT&MOVE
TECH

Catalogue

Product Functions	01
Product Diagram	01
Charging Method	02
Indicator Light Instructions	02
App Installation	03
Usages	04
Precautions	10
Product Parameters	11
Package Contents	11
Warranty	12


Product Functions

The Smart Grip Ball is used to do hand training and exercises with real-time data monitoring. Individuals who need hand rehabilitation, hypertension control, and improvements in muscular, cardiovascular and neurological health, such as those with concerns about Alzheimer's, Parkinson's, ALS, stroke, etc. can benefit from regularly using this device.

Main Functions:


1. Test grip strength;
2. Grip or wave to play games;

Product Diagram

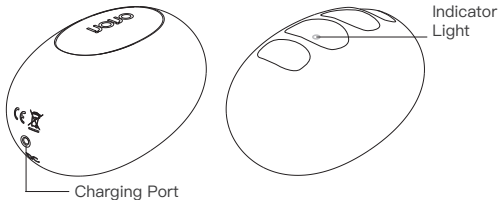

Charging Method

Plug the DC cable into the charging port to charge the battery. The grip ball vibrates, and the power is connected.

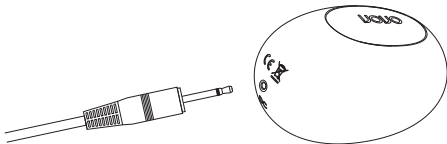

Indicator Light Instructions

	Light Effect	Status
Standby/Work	Flash	Power on
	Always on	Connection successful
Battery Indicator	Slow flash	Charging
	Always on	Fully charged
	Fast flash	Low battery

App Installation

Search “oriori Grip” in the Android Application Market, the iOS App Store, the HUAWEI Application Market to download and install the App.

App downloaded via other channels may be unusable. Before pairing with your equipment, make sure your equipment meets the following requirements:

- iPhone 5S and above with iOS 10.0 and above system;
- Android 5.0 and above with Bluetooth 4.0



App Store



Google Play



HUAWEI

Usages

1. Download and install “ oriori Grip ” App;
(Make sure the grip ball has enough battery power.)




2. Turn on the phone's Bluetooth;




ON / OFF




3. Open “oriOri Grip” App;
(Sign up when using for the first time.)

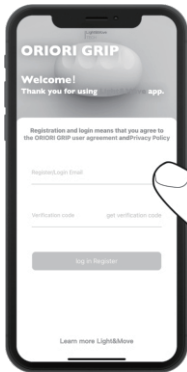

4. Bind your grip ball;
(Wave the grip ball until it vibrates.)

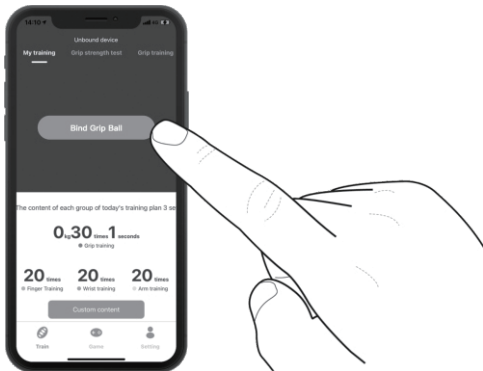

5. Find your grip ball serial number, select and connect;
(Click "Search again" if not found.)

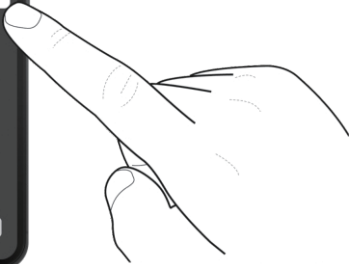

6. Confirm binding by grip.

(Make sure the grip ball has enough battery power.)



Precautions

- Only use accessories approved by the product manufacturer and compatible with this model. If you use other types of accessories, you may violate the warranty terms of the device and the applicable regulations of the country where the product is sold, and unexpected safety accidents may be caused.
- Do not expose the product and its accessories to high temperature or around heating devices, such as sunlight, heaters, microwave ovens, ovens or water heaters. Please use this product and its accessories under 0°C–45°C. If the environment temperature is too low or too high, product failure may occur.
- When charging is completed or this product is not charged, please disconnect the charger from the product and unplug the charger from the power socket.
- This product is equipped with a non-removable built-in battery. Please do not replace the battery yourself to avoid damage to the battery or the product. It can be replaced by an authorized service center.

- Please dispose of this product, batteries and other accessories in accordance with local regulations, and do not dispose of them as domestic waste.
- The elderly must use the grip ball in the company of others.
- If the grip ball is accidentally exposed to water, please wipe it and dry it before use. Soaking in the water or using it in the water is prohibited.

Product Parameters

Product Name: motion series

Size: 42×42×67.5 (mm)

Main Materials: Silicone

Rated Voltage: 3.7V

Rated Current: 5mA

Package Contents

Smart grip ball×1

Charging cable×1

Pouch×1

Instruction manual×1

Warranty

1. If any malfunction or quality problem is found within 1 year after the product is received, Please contact the distributor or customer service for after-sales support.
2. If it is found that breakage or damage is caused by force majeure factors such as improper use, man-made damage, natural disasters, such breakage or damage will not be covered by the warranty service.



For more information
please visit: <https://www.lightmovetech.com/>
Customer Service: info@lightmovetech.com